



MAIN FEATURES & BENEFITS

- Excellent material resistance towards aggressive media
- Inherently hydrophobic PTFE membrane
- Sterile grade at 0.2 µm (HIMA/ASTM)
- High flow rates
- Biologically inert
- Approved for Food Contact Use acc. to CFR Title 21 & 1935/2004/EC

PRODUCT DESCRIPTION

The LifeTec PT N filter element is a sterile grade, pleated high performance PTFE membrane filter. It provides the greatest assurance of filtration performance, stability and durability against chemicals even under severe process conditions.

The PTFE filter medium is inherently hydrophobic with a highly porous membrane structure. This ensures high flow rates and a high absorption of particles during the entire service life time. The end caps and the PTFE membrane are thermally welded without the use of binders. This results in an integral filter cartridge which provides maximum durability range against chemicals with minimal extractables.

This extremely durable design maintains consistent porosity and impurity retention throughout its service life without shedding or unloading contaminations.

All components meet the EU and USA requirements for Food Contact Use in accordance with CFR (Code of Federal Regulations) Title 21 and 1935/2004/EC. The filter element is manufactured in accordance with the manufacturing requirements, has no migration of filter media, is non-fibre releasing and is thermally welded.

INDUSTRIES



- Food



- Dairies



- Pharmaceutical



- Chemical

Donaldson Filtration Deutschland GmbH

Büssingstraße 1
42781 Haan • Germany
Tel. +49 2129 569 0
Fax +49 2129 569 100
CAP-de@donaldson.com
www.donaldson.com

APPLICATIONS

The sterile grade LifeTec PT N membrane filter is designed and developed for the following applications:

Purification of chemicals:

- Highly concentrated Acids
- Highly concentrated Bases
- Oxidizing reagents
- Complexing agents
- Alcohols, Aldehydes
- Etchants
- Chlorinated and fluorinated solvents
- Esters and Ketones
- Photolithographic Liquids

Filtration of air and gases:

- Compressed Air
- Fermentation Air
- Tank Ventilation
- Technical Gases

MATERIAL COMPLIANCE EU

The Donaldson LifeTec PT N filter element meets the guideline for Food Contact Use as given in European Regulation (EC) Number 1935/2004. All polymeric components (Polypropylene, PTFE) meet the requirements of EU Directive EC/10/2011 relating to plastic materials and articles intended to come into contact with foodstuffs.

Migration tests have been carried out in simulants after flushing or in flow conditions.

For specific details on the O-rings, please contact your Donaldson Sales Engineer.

RETENTION RATES (According to HIMA Challenge per ASTM)

Filter Grade	Microorganism	LRV / cm ²
LifeTec PT N 0.2 µm	Serratia Marcescens	> 7
	Pseudomonas diminuta	> 7

MATERIAL COMPLIANCE USA

All components of the LifeTec PT N filter element are FDA listed for food contact use in the Code of Federal Regulations (CFR), Title 21

Filter Materials		CFR Title 21
Membrane	PTFE	§ 177.1550
Upstream Support	Polypropylene	§ 177.1520
Downstream Support	Polypropylene	§ 177.1520
Outer Guard	Polypropylene	§ 177.1520
Core	Polypropylene	§ 177.1520
End Caps	Polypropylene	§ 177.1520
O-Rings	EPDM	§ 177.2600
	Silicone	§ 177.2600
Sealing Method	Thermal Bonding	

QUALITY TEST

All products have been inspected and released by Quality Assurance as having met the following requirements:

- All sterile filters are integrity tested to verify compliance with established quality and design specifications and to assure consistent and reliable performance.
- All LifeTec PT N filter elements are completely staged, assembled, tested and packaged in Class 7 clean room facility, whose Quality Management System is approved by an accredited registering body to the appropriate ISO 9001 Quality Systems Standard.


PRODUCT SPECIFICATIONS

Product Specifications				
Filter Grade	0.2 µm (Retention Rates LRV >= 7 cm ²)			
Filtration Surface	0.85 m ² per 250 mm element (10")			
Maximum Differential Pressure	Operating temperature		Differential pressure	
	°C	°F	bar	psi
	38	100	5.5	80
	66	150	4.1	60
	82	180	2.1	30
Cumulative Steaming Time*	121°C – 125°C (30 minutes) Saturated Steam (Forward Flow) up to 100 cycles			

* Figures are based on lab tests to evaluate steaming resistance. Filter elements need to be checked in actual use. Contact Donaldson for recommended Autoclaving/Steaming procedures.

For the filtration of aqueous solutions the LifeTec PT N membrane filter has to be pre-wetted with a suitable liquid of low surface tension (e.g. IPA).


FLOW CHARACTERISTICS

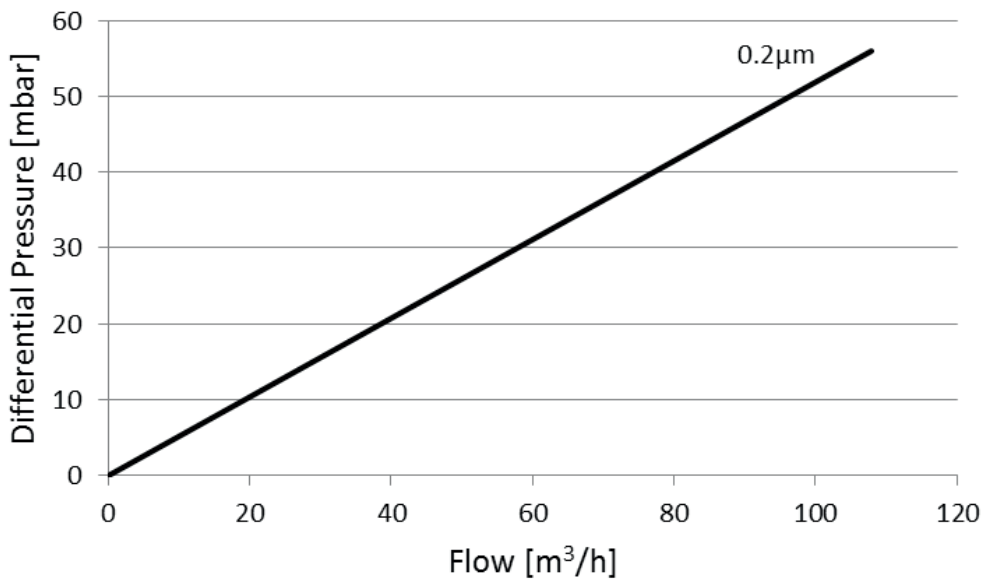

LifeTec PT N
10", air, 25°C, 1 bar absolute


INTEGRITY TESTING

Water Intrusion Test	
Filter Grade	Maximum Diffusion Value
LifeTec PT N 0,2 µm	0.60 ml/min per 10" element after 10 min @ 2.0 bar (29 psi)

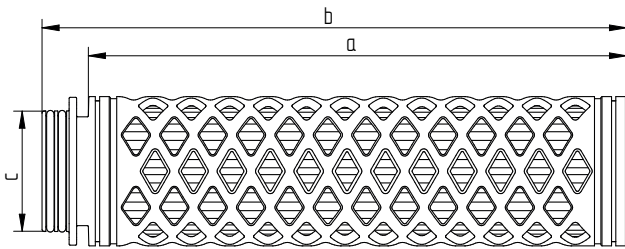


Dimensions (CODE 2 connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	253	10.0	274	10.8	56	2.2
20"	495	19.5	516	20.3	56	2.2
30"	737	29.0	758	29.8	56	2.2
40"	979	38.5	1000	39.4	56	2.2

CODE 2: 2 x 226 o-rings, bayonet 2 locking tabs, flat end cap, integrated reinforcement ring

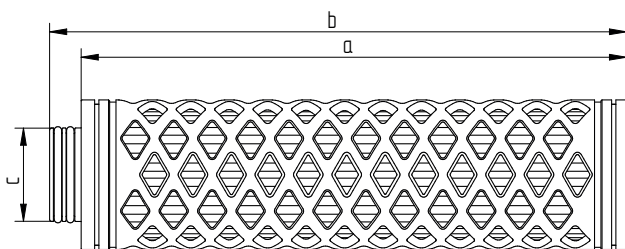


Dimensions (CODE 3 connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	256	10.1	271	10.7	44	1.7
20"	498	19.6	513	20.2	44	1.7
30"	740	29.1	755	29.7	44	1.7
40"	982	38.7	997	39.3	44	1.7

CODE 3: 2 x 222 o-rings, plug connection, flat end cap, integrated reinforcement ring

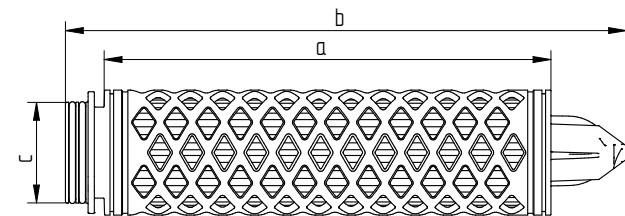


Dimensions (CODE 7 connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	251	9.9	315	12.4	56	2.2
20"	493	19.4	557	21.9	56	2.2
30"	735	28.9	799	31.5	56	2.2
40"	977	38.5	1041	41.0	56	2.2

CODE 7: 2 x 226 o-rings, bayonet 2 locking tabs, locating fin, integrated reinforcement ring

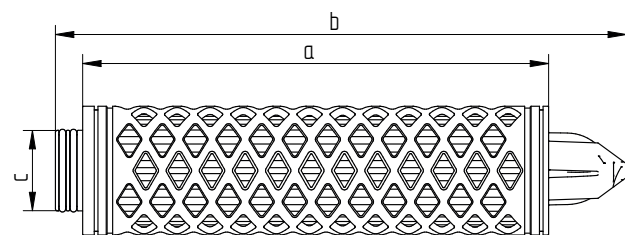


Dimensions (CODE 8 connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	254	10.0	311	12.2	44	1.7
20"	496	19.5	553	21.8	44	1.7
30"	738	29.1	795	31.3	44	1.7
40"	980	38.6	1037	40.8	44	1.7

CODE 8: 2 x 222 o-rings, plug connection, locating fin, integrated reinforcement ring

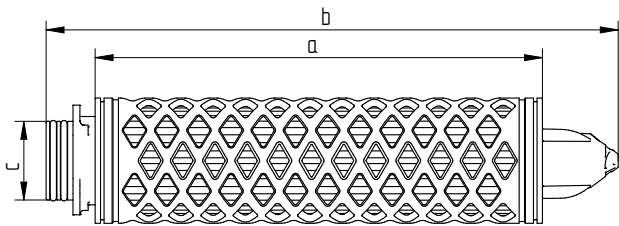


Dimensions (CODE 9 connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	250	9.8	320	12.6	44	1.7
20"	492	19.4	562	22.1	44	1.7
30"	734	28.9	804	31.7	44	1.7
40"	976	38.4	1046	41.2	44	1.7

CODE 9: 2 x 222 o-rings, bayonet 3 locking tabs, locating fin, integrated reinforcement ring

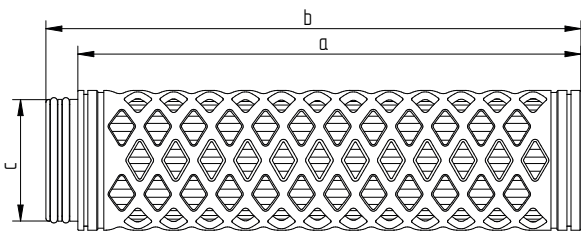

Dimensions (UF connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	252	9.9	268	10.6	61	2.4
20"	494	19.4	510	20.1	61	2.4
30"	736	29.0	752	29.6	61	2.4

CODE UF: 2 x 226 o-rings, plug connection, flat end cap, integrated reinforcement ring

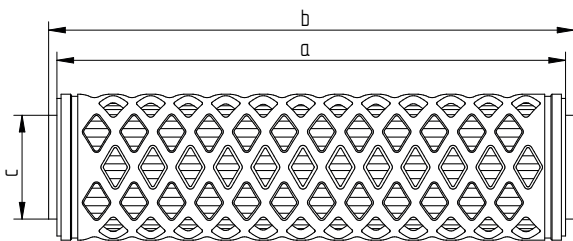

Dimensions (DOE connection)						
Size	a		b		c	
	mm	inch	mm	inch	mm	inch
10"	244	9.6	250	9.8	50	2.0
20"	500	19.7	506	19.9	50	2.0
30"	754	29.7	760	29.9	50	2.0
40"	1008	39.7	1014	39.9	50	2.0

DOE: Double open end with EPDM gaskets

Other end cap configurations on request.

- Integrity test to be done by Water Intrusion Test.
- For information on test equipment or test services, please contact your Donaldson Sales Engineer and visit our website at www.donaldson.com!